UltraSeal XT[®] *plus*[™]

Pit & Fissure Sealant
and Flowable Composite

Description:

- Thixotropic, resin sealant that penetrates into occlusal grooves, when used with the Inspiral[®] Brush tip that allows for shear thinning effect of the sealant.
- 58% filled for optimum strength.
- UltraSeal XT plus in conjunction with PrimaDry[®] drying and priming agent virtually eliminates microleakage.
- Radiopaque.

Indications:


Use UltraSeal XT plus for prophylactic sealing of pits and fissures.


Thixotropic, 58% filled UltraSeal XT plus flows out of a syringe through the helical channel of the Inspiral Brush tip.

Preliminary Procedure:


1. Remove luer lock cap from Ultra-Etch[®] 35% phosphoric acid syringe. Firmly twist a disposable Inspiral Brush tip or Blue Micro[®] tip onto syringe.
2. Remove luer lock cap from 1.2 ml PrimaDry syringe. Firmly twist either a Black Mini[®] Brush or a Blue Micro tip onto syringe.
3. Remove luer cap from UltraSeal XT plus syringe. Firmly twist an Inspiral Brush tip or a Black Micro[®] FX tip onto syringe.
4. Clean fissures with a prophy-jet type device, high pressure air abrasion or a thin diamond bur. We recommend the Fissure Prep Bur to remove enamel pellicle, etc.; If using a prophy-jet type device with sodium bicarbonate, it is advised to etch twice. The first application neutralizes the sodium bicarbonate; the second conditions the enamel.
5. Hold syringe in palm of hand for optimum control rather than in conventional manner with thumb on plunger.

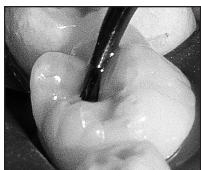
Etch

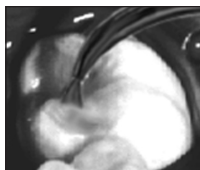Dry and Prime

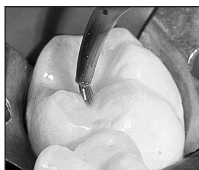Seal

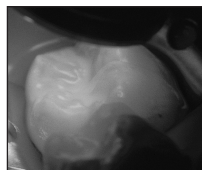Cure

Procedure:

Pit & Fissure Sealants

Etch:

1. Isolate teeth with rubber dam or appropriate absorbent devices such as cotton rolls, etc. Moisture contamination either in liquid or vapor form just prior to resin placement will compromise bond and seal.¹
2. Deliver Ultra-Etch to fissure using the Inspiral Brush tip or the Blue Micro tip, according to preference. Etch with Ultra-Etch for 15-20 seconds, then rinse with full air/water spray and dry.

Dry and Prime:

Apply PrimaDry and leave 5 seconds. Dry by gently blowing area with moisture-free and oil-free air (don't rinse).²

Note: PrimaDry significantly enhances bond and seal of UltraSeal XT plus. Remember, fissures are not empty spaces. Organic matter in the fissures can act like a sponge, supplying moisture contamination. PrimaDry rapidly volatilizes moisture content of fissures.

Note: Most dental 3-way syringes deliver moist contaminants with air. Use Aqua✓® to evaluate moisture content of air in 3-way syringe.

Note: If applying to dentin, use a dentin bonding agent such as PQ Clear in place of PrimaDry for both dentin and enamel. Do not use PrimaDry as a dentin bonding agent.

Seal:

The Inspiral bristle length can be customized by grasping bristles between thumb and finger and moving them in or out. A shorter length is desirable to “stiffen” the brush; We recommend about

1.2mm in length. Express a small drop of sealant at the end of the Inspiral Brush tip. Using a scrubbing motion, apply to the deepest anatomy of the fissure. The resin is then worked into the fissures with the tip of the fibers. The Black Micro FX tip is an alternative delivery tip. It is preferable to not place too much resin.

Light Cure:

Cure 20 seconds if using a standard curing light, 10 seconds if using a high energy light. (Use a quality Halogen or the quality Ultra-Lume® LED lights.) Adjust occlusion.

Precautions:

1. Replace luer lock caps when syringes are not in use or use TipSoc™ covers.
2. All 1.2 ml tips are disposable: use one per patient.
3. Clean and disinfect 1.2 ml syringes between patients.
4. Never attempt to deliver any Ultra-Etch or sealant through the Black Mini Brush tip.
5. Always twist tips securely onto syringes. Use slow, steady pressure when expressing materials through brush tips. Always check the flow of each material on a mixing pad before intra-oral use. Never force syringe plungers as this could result in the tip blowing off and an uncontrolled discharge of material. If there is resistance to flow, replace with a fresh tip and test again.
6. Remember that any strong acid such as phosphoric acid can burn and irritate tissues. Be careful not to splash into eyes, onto mucous membranes, skin or any other area. Confine to areas being etched only. Any other area of exposure should be immediately rinsed thoroughly with water.
7. PrimaDry is flammable. Handle with care.
8. Following evaporation of PrimaDry, a slight residue will remain. This will re-dissolve should more PrimaDry be applied. To prevent drying of volatile materials such as PrimaDry, a TipSoc can be applied.
9. Any resin can be sensitizing to the patient. Do not use on patients with known sensitivity/allergy to resins. Wash exposed area completely to remove material. If dermatitis, rash or any other allergic reaction is observed, consult a physician.

¹ Lee, J., "Microleakage of New Light-Cured Sealants," Practical Periodontics and Aesthetic Dentistry, Supplement No. 1, pp. 19-23, June/July 1995.

² Gene L. Dickinson, D.D.S., M.S., "Microleakage Evaluation of UltraSeal XT Pit and Fissure Sealant with PrimaDry Drying Agent," Medical College of Georgia, December 8, 1992. Available upon request.

³ Studies (two examples below) and clinical experience show that quality of light varies from manufacturer to manufacturer.

Important: Effectiveness of light and light conducting fibers, tubes, etc., decreases as unit ages. Light cure units should be checked to be sure that maximum light is being delivered for proper curing. Lower intensity lights may require additional curing time. Check with a curing radiometer.

A. Dugan, William T. and Hartleb, John H. "Influence of a Glutaraldehyde Disinfecting Solution On Curing Light Effectiveness." General Dentistry 37: 40-43, 1989.

B. Clinical Research Associates Newsletter, 12: 3, February, 1988.

UltraSeal XT[®] plus[™]

Fissure sealant and flowable composite

DE-Fissurenversiegler & fließfähiges Composite
FR-Scellement de fissure et composite fluide
NL-Fissuursealend vloeibaar composietmateriaal
IT-Sigillante per solchi e composito scorrevole
ES-Sellador de fisuras y composite fluido
PT-Selante de fissuras e composto dispersável
SV-Spricktätningsmedel med flytande sammansättning
DA-Fissurforsøgning og komposit, der kan flyde
FI-Fissuuroiden pinnoite ja valuva yhdistelmämuovi
EL-Υλικό έμφραξης σχισμών και υλικό αποκαταστάσεων υψηλής ρευστότητας

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



Inspiral[®] Brush Tip

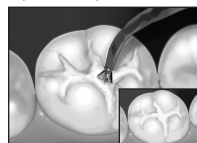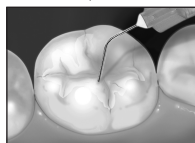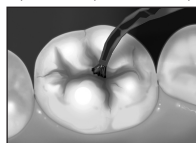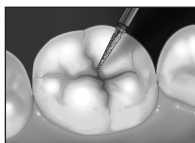EN Freshen enamel using Fissure Prep Bur, air abrasion, etc.

Ultra-Etch[®] 15 seconds.
Air/water wash 5 seconds. Air dry.

PrimaDry[®] 5 seconds; gently air dry. Significantly enhances seal!

Work UltraSeal XT plus into surface. Light cure 20 seconds.

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EN Precautions:

1. For Professional use only.
2. Review instructions, precautions, and MSDS before beginning treatment. Use only as directed.
3. Keep products out of heat/sunlight.
4. Avoid skin exposure to resins.
5. Isolate strong chemicals to area of treatment.
6. Confirm that patient has no known allergies to treatment materials.
7. Test flow of materials from syringe and tip before using intraorally.
8. Never force syringe plungers.

9. Use Ultradent® Syringe Covers and/or clean and disinfect syringes between patients.
10. Use standard curing lights with all light polymerized/activated materials.
11. All tips are disposable.
12. Never attempt to deliver Ultra-Etch or UltraSeal XT plus through the Black Mini Brush tip.
13. PrimaDry is flammable. Handle with care.
14. To prevent drying of volatile materials such as PrimaDry, replace luer lock caps when syringes are not in use or use TipSoc covers.

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ULTRADENT PRODUCTS, INC. **MATERIAL SAFETY DATA SHEET**

Controlled Product Yes ☐ No ☒

SECTION I - IDENTIFICATION OF THE SUBSTANCE/ PREPARATION AND OF THE COMPANY/ UNDERTAKING

Material Name UltraSeal XT [®] plus [™]	Manufactured by Ultradent Products, Inc. 505 West 10200 South South Jordan, Utah 84095 USA
Material Use Pit and Fissure Sealant and Flowable Composite.	Emergency Contact 1-800-552-5512 or 1-801-572-4200 (Outside USA)
Date August 29, 1997	EC Representative
Prepared by: Steve Jensen, Chemist	

SECTION II - COMPOSITION/ INFORMATION ON INGREDIENTS

Active Ingredients	Approximate Concentration %	LD ₅₀ /LC ₅₀ Specify Species and Route	C.A.S. N.A. or U.N. Numbers
Bis-GMA	proprietary	Not established	1565-94-2

SECTION III - HAZARDS IDENTIFICATION

Adverse Human Effects and Symptoms May cause sensitization.

SECTION IV - FIRST AID MEASURES

Inhalation N/A	Eye Contact Flush with lots of water. Call a physician if needed.
Skin Contact Wash with soap and water.	Ingestion In large amounts, induce vomiting.

SECTION V - FIRE-FIGHTING MEASURES

Flammability If Yes, Under Which Conditions: Avoid contact with open flame.
Yes ☒ No ☐

Means of Extinction CO₂, dry chemical.

Protective Equipment for Firefighters Standard.

SECTION VI - ACCIDENTAL RELEASE MEASURES

Personal & Environmental Precautions None.

SECTION VII - HANDLING AND STORAGE

Methods for Cleaning Wipe up as much as possible with dry cloth.
Clean remainder with alcohol.

Handling Precautions None.

Storage Requirements Refrigerate when not in use.

SECTION VIII - EXPOSURE CONTROLS/ PERSONAL PROTECTION

Respiratory Protection N/A	Eye Protection Protective eyewear.
Hand Protection Latex gloves.	Skin Protection Lab coat.

SECTION IX - PHYSICAL & CHEMICAL PROPERTIES

Appearance Tinted liquid.	Solubility in Water Insoluble.	Melting Point/Range N/A
Odor Slight acrylic odor.	pH N/A	Other Data None.

SECTION X - STABILITY AND REACTIVITY

Chemical Stability Conditions and Materials to Avoid Light exposure will initiate polymerization.

Hazardous Decomposition Products CO₂

SECTION XI - TOXICOLOGICAL INFORMATION

Route of Entry ☒ Skin Contact ☒ Skin Absorption ☒ Eye Contact ☐ Inhalation Acute ☐ Inhalation Chronic ☐ Ingestion
Effects and Symptoms Skin: Slight irritation. Eye: Slight irritation.
Chronic Effects and Symptoms of Long-Term Exposure Skin: Low possibility of sensitization upon prolonged exposure. Eye: Slight irritation.
Sensitization Slight possibility.

Irritation of Product Slight.

SECTION XII - ECOLOGICAL INFORMATION

Biodegradability N/A

SECTION XIII - DISPOSAL CONSIDERATIONS

Method of Waste Disposal

Further Details None.

In accordance with government regulations.

SECTION XIV - TRANSPORT INFORMATION

Not regulated.

SECTION XV - REGULATORY INFORMATION

Not defined.

SECTION XVI - OTHER INFORMATION

FOR DENTAL USE ONLY. Use as directed.

The information and recommendations are taken from sources (raw material MSDS(s) and manufacturer's knowledge) believed to be accurate; however, Ultradent Products, Inc., makes no warranty with respect to the accuracy of the information or the suitability of the recommendation and assumes no liability to any user thereof. Each user should review these recommendations in the specific context of the intended use and determine whether they are appropriate.

ULTRADENT PRODUCTS, INC. **MATERIAL SAFETY DATA SHEET**

Controlled Product Yes ☐ No ☒

SECTION I - IDENTIFICATION OF THE SUBSTANCE/PREPARATION AND OF THE COMPANY/UNDERTAKING

Material Name	Ultra-Etch®	Manufactured by	Ultradent Products, Inc. 505 West 10200 South South Jordan, Utah 84095 USA
Material Use	Etching of tooth enamel and dentin.	Emergency Contact	1-800-552-5512 or 1-801-572-4200 (Outside USA)
Date	June 10, 1997	EC Representative	
Prepared by:	Steve Jensen, Chemist		

SECTION II - COMPOSITION/INFORMATION ON INGREDIENTS

Active Ingredients	Approximate Concentration %	LD ₅₀ /LC ₅₀ Specify Species and Route	C.A.S. N.A. or U.N. Numbers
Phosphoric Acid	35%	LD ₅₀ 1530 mg/kg Oral – Rat	7664-38-2

SECTION III - HAZARDS IDENTIFICATION

Adverse Human Effects and Symptoms Irritates eyes and skin.

SECTION IV - FIRST AID MEASURES

Inhalation	N/A	Eye Contact	Flush immediately with water. Get medical attention.
Skin Contact	Flush with water.	Ingestion	Do not induce vomiting. Give large amounts of water or milk. Take an antacid. Call a physician.

SECTION V - FIRE-FIGHTING MEASURES

Flammability If Yes, Under
Yes ☐ No ☒ Which Conditions:

Means of Extinction N/A

Protective Equipment for Firefighters N/A

SECTION VI - ACCIDENTAL RELEASE MEASURES

Personal & Environmental Avoid contact with skin, eyes, and clothing. Prevent contamination of soil, drains, etc.

SECTION VII - HANDLING AND STORAGE

Methods for Cleaning Clean up with wet reg. Neutralize with baking soda.

Handling Precautions Use gloves and protective eyewear.

Storage Requirements Store at room temperature.

SECTION VIII - EXPOSURE CONTROLS/PERSONAL PROTECTION

Respiratory Protection N/A

Eye Protection Protective eyewear.

Hand Protection Latex gloves.

Skin Protection Lab coat.

SECTION IX - PHYSICAL & CHEMICAL PROPERTIES

Appearance Blue gel (Ultra-Etch); Light blue gel (Ultra-Etch AB)

Solubility in Water Partial. Melting Point/Range Not defined.

Odor Odorless. pH < 2

Other Data

SECTION X - STABILITY AND REACTIVITY

Chemical Stability Conditions and
Yes ☒ No ☐ Materials to Avoid

Hazardous Decomposition Products Oxides of phosphorus.

SECTION XI - TOXICOLOGICAL INFORMATION

Route of Entry ☒ Skin Contact ☐ Skin Absorption ☒ Eye Contact ☐ Inhalation Acute ☐ Inhalation Chronic ☒ Ingestion

Effects and Symptoms of Short-Term Exposure Skin: May cause irritation. Eyes: Severe irritation. Ingestion: Possible irritation.

Chronic Effects and Symptoms of Long-Term Exposure Skin: Mild irritation. Eyes: Damage (chronic exposure must be avoided). Ingestion: Irritation

Sensitization N/A

Irritation of Product Mild.

SECTION XII - ECOLOGICAL INFORMATION

Biodegradability Not defined.

SECTION XIII - DISPOSAL CONSIDERATIONS

Method of Waste Disposal In accordance with governmental regulations.

Further Details N/A

SECTION XIV - TRANSPORT INFORMATION

No limitations.

SECTION XV - REGULATORY INFORMATION

Not defined.

SECTION XVI - OTHER INFORMATION

FOR DENTAL USE ONLY. Use as directed.

The information and recommendations are taken from sources (raw material MSDS(s) and manufacturer's knowledge) believed to be accurate; however, Ultradent Products, Inc., makes no warranty with respect to the accuracy of the information or the suitability of the recommendation and assumes no liability to any user thereof. Each user should review these recommendations in the specific context of the intended use and determine whether they are appropriate.

ULTRADENT PRODUCTS, INC. **MATERIAL SAFETY DATA SHEET**

Controlled Product Yes ☐ No ☒

SECTION I - IDENTIFICATION OF THE SUBSTANCE/PREPARATION AND OF THE COMPANY/UNDERTAKING

Material Name	PrimaDry®	Manufactured by	Ultradent Products, Inc. 505 West 10200 South South Jordan, Utah 84095 USA
Material Use	Drying and Priming Agent.	Emergency Contact	1-800-552-5512 or 1-801-572-4200 (Outside USA)
Date	September 2, 1997	EC Representative	
Prepared by:	Steve Jensen, Chemist		

SECTION II - COMPOSITION/INFORMATION ON INGREDIENTS

Active Ingredients	Approximate Concentration %	LD ₅₀ /LC ₅₀ Specify Species and Route	C.A.S. N.A. or U.N. Numbers
Ethyl Alcohol	99%	LD ₅₀ 7500mg/kg Oral - Mouse LC ₅₀ 39g/m 3/4H Inhalation - Mouse (no toxic effect noted)	64-17-5

SECTION III - HAZARDS IDENTIFICATION

Adverse Human Effects and Symptoms Ethanol ingested in large volumes can be toxic.

SECTION IV - FIRST AID MEASURES

Inhalation	N/A	Eye Contact	Flush with lots of water.
Skin Contact	Wash with soap and water.	Ingestion	Give large amounts of water or milk.

SECTION V - FIRE-FIGHTING MEASURES

Flammability If Yes, Under Which Conditions: Exposure to spark, flame, or heat.
Yes ☒ No ☐

Means of Extinction Dry chemical, CO₂, water spray.

Protective Equipment for Firefighters Standard protection.

SECTION VI - ACCIDENTAL RELEASE MEASURES

Personal & Environmental Precautions Not defined.

SECTION VII - HANDLING AND STORAGE

Methods for Cleaning Wipe up with paper towel or cloth.

Handling Precautions Not defined.

Storage Requirements Store at room temperature. Keep out of heat and/or direct sunlight.

SECTION VIII - EXPOSURE CONTROLS/PERSONAL PROTECTION

Respiratory Protection	N/A	Eye Protection	Protective eyewear.
Hand Protection	Latex gloves.	Skin Protection	Lab coat.

SECTION IX - PHYSICAL & CHEMICAL PROPERTIES

Appearance	Clear liquid.	Solubility in Water	Completely.	Melting Point/Range	N/A
Odor	Ethyl alcohol odor.	pH	N/A	Other Data	Rapid evaporation rate. 99% volatile.

SECTION X - STABILITY AND REACTIVITY

Chemical Stability Conditions and Materials to Avoid Strong oxidizing agents.

Hazardous Decomposition Products CO₂

SECTION XI - TOXICOLOGICAL INFORMATION

Route of Entry	<input type="checkbox"/> Skin Contact <input type="checkbox"/> Skin Absorption <input checked="" type="checkbox"/> Eye Contact <input type="checkbox"/> Inhalation Acute <input type="checkbox"/> Inhalation Chronic <input checked="" type="checkbox"/> Ingestion
Effects and Symptoms of Short-Term Exposure	Eyes: Irritation.
Chronic Effects and Symptoms of Long-Term Exposure	Eyes: Irritation. Ingestion: Irritation (burning of throat).
Sensitization	N/A

SECTION XII - ECOLOGICAL INFORMATION

Biodegradability Yes.

SECTION XIII - DISPOSAL CONSIDERATIONS

Method of Waste Disposal In accordance with government regulations.

Further Details None.

SECTION XIV - TRANSPORT INFORMATION

Not regulated.

SECTION XV - REGULATORY INFORMATION


Not defined.

SECTION XVI - OTHER INFORMATION

FOR DENTAL USE ONLY. Use as directed.


The information and recommendations are taken from sources (raw material MSDS(s) and manufacturer's knowledge) believed to be accurate; however, Ultradent Products, Inc., makes no warranty with respect to the accuracy of the information or the suitability of the recommendation and assumes no liability to any user thereof. Each user should review these recommendations in the specific context of the intended use and determine whether they are appropriate.

UltraSeal XT plus

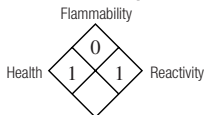
MSDS Included.
Shelf Life: 24 months

Ultra-Etch

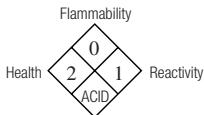MSDS Included.
Shelf Life: 48 months

PrimaDry

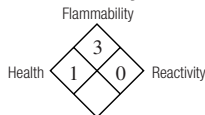MSDS Included.
Shelf Life: 24 months

HAZARD RATING

- 4 = Severe
- 3 = Serious
- 2 = Moderate
- 1 = Slight
- 0 = Minimal

Ultradent syringes have an expiration date stamped on the side of the syringe consisting of one letter and three numbers. The letter is a lot number used for manufacturing purposes and the three numbers are the expiration date. The first two numbers are the month, and the third number is the last number of the year.

U.S. federal law restricts the sale of this device by or on the order of a dentist.
Keep out of reach of children.
Refrigerate if not used on a daily basis.

For immediate reorder and/or complete descriptions of Ultradent's product line, refer to Ultradent's catalog or call Toll Free **1-800-552-5512**.
Outside U.S. call (801) 572-4200 or visit www.ultradent.com.

© Copyright 2004 Ultradent Products, Inc. All Rights Reserved.
U.S. Patent Nos.: 4,997,371; 5,246,371; 5,269,684; and 5,425,641. Other U.S. and International Patents Granted and Pending.

ULTRADENT
PRODUCTS, INC.

Manufactured by: Ultradent Products, Inc.
505 West 10200 South, South Jordan, Utah 84095, USA

#10907.10 091504